**ZAMST EK-1 / EK-3 / EK-5 (KNEE SUPPORT) INSTRUCTIONS FOR USE**

For safety reasons, be sure to read this user manual before use. Also, store it in a location where you can refer to it at any time.

**Read this first:**

- This product can be used for the left or right side.
- Depending on usage, this product may not be able to provide support.
- This product is not intended to be used for medical treatment.

**When using the product:**

- Creases that may seem to be a flaw may appear in fabric. However, this does not affect product quality.
- Due to the characteristics of the materials used in this product, the color may change under certain storage conditions.
- The product should be used with care as hook-and-loop fasteners may cause clothing to become frayed or damaged.
- When storing this product, do not flatten the product by applying excessive weight to it. Otherwise, the hook-and-loop fastener may not work effectively.

**Precautions for use - Be sure to read this:**

Safety precautions are divided into the following sections, according to the level of danger involved.

**⚠ WARNING:** Indicates a potentially hazardous situation, which may result in serious personal injury if not avoided.

**⚠ CAUTION:** Indicates a potentially hazardous situation, which may result in minor personal injuries and/or property damage if not avoided.

**⚠ WARNING**

- To prevent abnormalities or worsening of symptoms, people who suffer from the following conditions should refrain from using the product:
  - Those predisposed to allergies or with sensitive skin.
  - Those who have a fracture, cut, numbness, swelling, eczema, rash or other abnormalities in the area where the product will be worn.
- If the product is not used correctly, it may cause circulatory problems or injury. Be sure to follow the warnings below and use the product correctly:
  - The product should only be used on knees.
  - Be sure to use the correct size product.
  - Do not use the product while sleeping.
  - Follow the attachment steps in order to use this product correctly.
  - Do not wrap the main body of the product or the straps too tightly.
  - If the product loosens or comes off during use, be sure to refasten it. If wearing the product for a long period of time, take extra care to maintain appropriate compression force, and readjust the tightness as necessary.
  - If the product shows any signs of wear or damage, or the hook-and-loop fastener no longer works well, do not use the product.
  - Do not attempt to modify or repair the product.
- If any abnormalities, such as an injury, fracture, pain, numbness, swelling, eczema, etc. occur when using this product, immediately stop using it, and consult a physician. Continued use may worsen symptoms.

**⚠ CAUTION**

- Hard materials are used for the hook-and-loop fastener. In certain rare instances, fingers or skin may be damaged, so use the product carefully.
- Sweat and friction may cause the color of this product to fade or be transferred to other materials.
- Make sure the hook-and-loop fastener does not attach to clothing. If it does, clothing may become frayed or be damaged.

**Care and cleaning**
**⚠ CAUTION**


- Do not wash this product with light colored items. By doing so, color may fade or be transferred to other materials.
- If washing the product together with clothing, etc., be sure to use a laundry net in order to prevent the clothing from becoming frayed or be damaged.

- This product can be washed in a washing machine. After washing, hang it to dry in the shade.
- When washing this product, attach the hook-and-loop fastener to the main body. Attempting to wash this product without closing the hook-and-loop fastener may weaken its fastening capability.
- Do not use an iron or chlorine bleach, however, oxygen bleach can be used.
- Do not tumble dry. Doing so may damage the product.
- This product cannot be dry cleaned. Wet cleaning is allowed, however, avoid tumble drying.


# How to apply

- Before attaching this product, carefully read through the "Precautions for use - Be sure to read this" section. Then follow the steps provided to ensure the product is correctly attached.
- The steps below use the right leg for reference.
- These steps show how to apply EK-3. Use the same steps to apply EK-1 or EK-5.


1. Place the center of the product at the back of your knee with the correct side up. Straighten your knee and fasten the upper hook-and-loop fastener B above your kneecap.

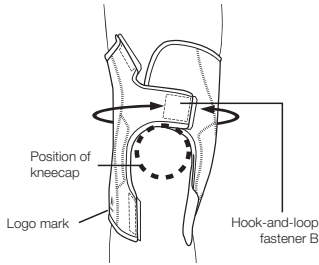
2. Fasten hook-and-loop fastener C below your kneecap.

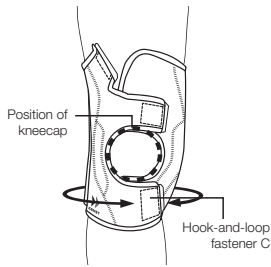
3. Fasten hook-and-loop fastener A.

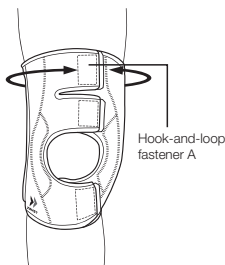4. Adjust hook-and-loop fasteners A, B, and C so that you will be able to maintain an appropriate compression force.

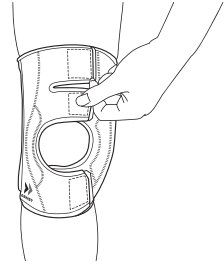
5. For EK-5 Only: Fasten the femoral cross straps at the back of your thigh.

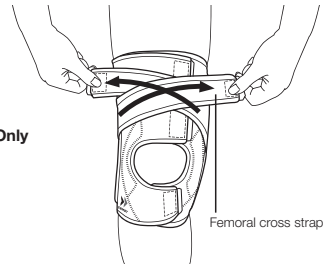**For EK-5 Only**

6. Pull the product up and down to eliminate slack.

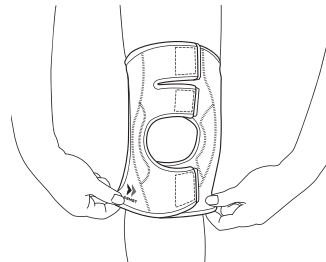7. Make sure that the product is not loose or incorrectly attached. If you feel any discomfort, repeat the steps from the beginning to reattach the product.

## Disclaimer:

- \* The product specifications and appearance are subject to change without notice for the purpose of improvement.
- \* When disposing of this product, be sure to follow all disposal ordinances and regulations as specified by the local government.
- \* We take all possible measures to ensure product quality, however, if you find any defects, etc. in our products, please contact Nippon Sigmex Co., Ltd.

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